



THIS
IS A
babyLance™
BABY

babyLance™ Blade and Incision Dimensions Guide

Designed for neonatal nurses...with neonatal nurses!

babyLance™ safety heelsticks have earned worldwide popularity for their **smooth and accurate incisions, ease of use and safety.**



Americas
3883 Rogers Bridge Road NW
Suite 501
Duluth, GA 30097
Tel: +1 770 448 9493

Asia
15 Hoe Chiang Road #12-02
Tower Fifteen
Singapore 089316
Tel: +65 63451588


Europe
Tel: +44 330 122 2630





MediPurpose™
medipurpose.com/babylance

Incision Dimensions¹

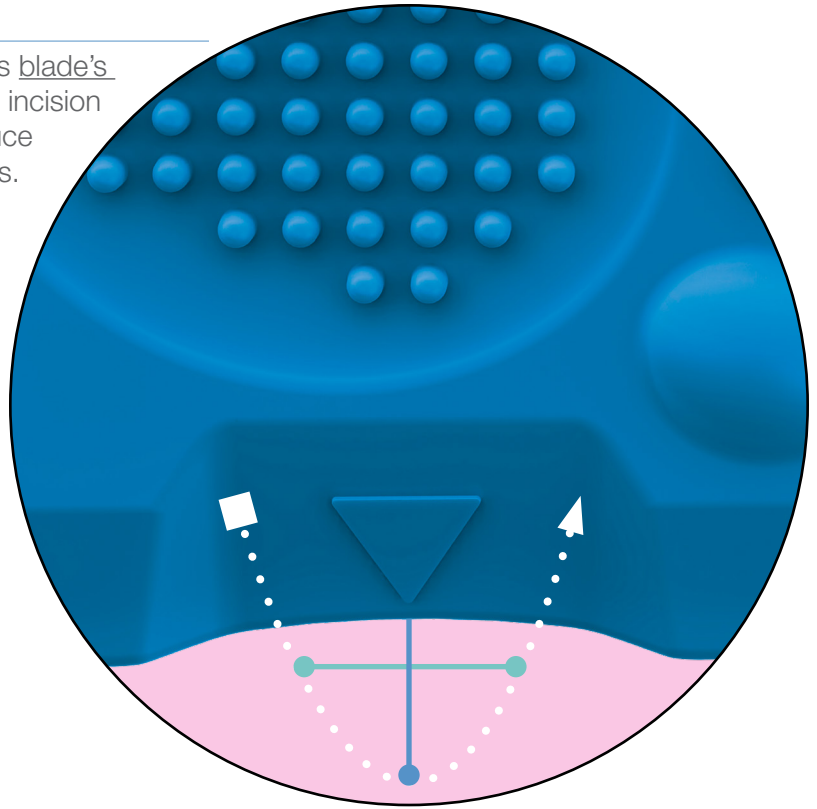babyLance's proprietary internal spring activates its blade's swift pendulum action—delivering an ultra-smooth incision that complies with CLSI LA4-A5 guidelines to reduce trauma to an infant's delicate subcutaneous tissues.

The **gently arcing incision** (white dotted line) is measured by the **length** (green horizontal line) and the **maximum depth** (blue vertical line).




babylance [™] Product	Incision Length ²	Incision Depth ³
 Micropreemie (BLM)	2.2 mm	0.60 mm
 Preemie (BLP)	3.0 mm	0.85 mm
 Newborn (BLN)	3.0 mm	1.0 mm
 Toddler (BLT)	4.0 mm	1.4 mm

Blade Dimensions¹

Blade Width: 3.0 mm

Blade Thickness: 0.15 mm

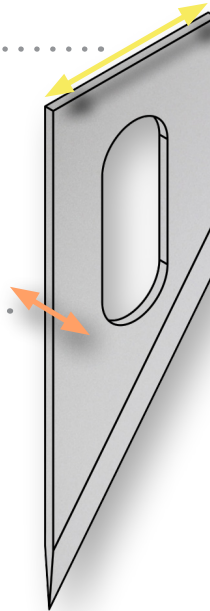

¹ Blade illustration and incision profile scales are enlarged for demonstration purposes only.

² +/- 0.5 mm

³ +/- 0.2 mm

Americas

3883 Rogers Bridge Road NW
Suite 501
Duluth, GA 30097
Tel: +1 770 448 9493

Asia

15 Hoe Chiang Road #12-02
Tower Fifteen
Singapore 089316
Tel: +65 63451588

Europe

Tel: +44 330 122 2630

MediPurpose[™]

medipurpose.com/babylance